


Firenze Bedhead Panel

www.lipmedical.com

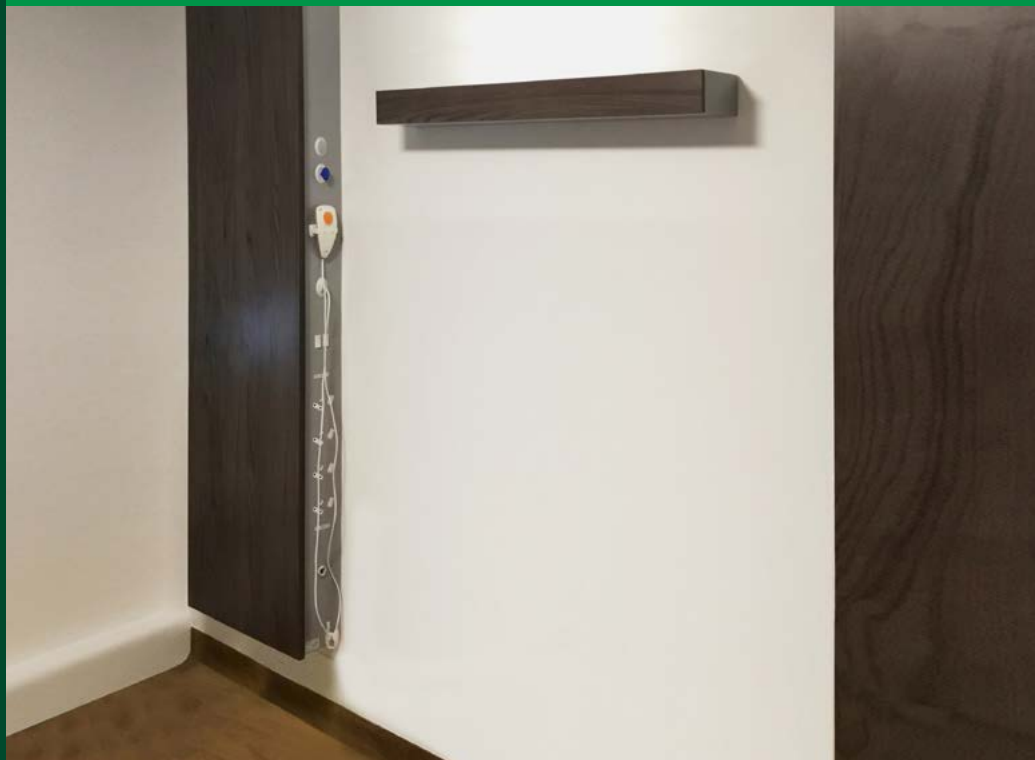

The future looks good.

Specifications:

Firenze Bedhead Panel

The Firenze offers a stylish compact solution to supply services to the patient bed space. The unit can be supplied with extruded aluminium channels down either one or both sides to provide electrical services and data with the medical gas outlets mounted in the HPL laminate front fascia. The standard width of the unit is 500mm and length 2000mm but these can be tailored to suit individual ceiling heights and client requirements. The front fascia panel is manufactured from high pressure laminate which can be selected from a standard range of finishes that includes block colours or woodgrains. Patterns or printed high resolution image graphic panels can also be provided to client choice and interior design. The Firenze bedhead units are provided as fully assembled units and can be provided either fully pre-wired and pre-gassed or to be gassed and wired on site.

Side Channels Can Be Fitted With :

- Single 13A socket outlets (switched or Isolated earth unswitched type)
- Equipotential Earth bonding jacks
- RJ45 Data Outlets
- Headphone Socket
- Nurse Call Panel (provision for any nurse call manufacturer provided)
- Lighting switches (to suit lighting design)
- USB Charging outlets

Fascia:

Up to 6 Medical gas outlets, Pneumatech East or Medaes GEM Shield as standard

Lighting:

Low level lighting with daylight sensors can be fitted to the underside of the units to provide ambient lighting during the evening

Dimensions:

Standard sizes - Width 500mm, Depth 155mm, Length 2000mm

T: 0113 251 5119

F: 0113 251 5100

enquiries@lipmedical.com

Offices

1200 Century Way, Thorpe Park Business Park, Colton, Leeds, LS15 8ZA
Unit 60 AML Zone 2, Deeside Industrial Estate, Third Avenue, Deeside, Flintshire, CH5 2LA

