

Mediline Horizontal Bedhead Unit

www.lipmedical.com



MANUFACTURED IN THE UK




LIP MEDICAL
HEALTHCARE DESIGN

Mediline Horizontal Bedhead Unit

The Mediline horizontal bedhead unit offers an ergonomic, compact stylish solution to deliver essential services to the patient. The units are configured to suit each individual installation or project requirement with modules supplied fully assembled as standard. Electrical sockets, nurse call, data and medical gases are fitted to meet each individual client specification.




Electrical sockets are fitted to individual accessory plates and can be fully wired and tested with essential/non essential sockets available for standard ward areas and blue isolated earth sockets for A&E, critical care or operating departments.

Medical gas terminal units from all manufacturers can be accommodated with modules available as either fully pregassed and tested systems or alternatively, mounting plates can be included for installation of the gas outlets by the medical gas contractor on site.

Any type of nurse call can be easily integrated with provision for a bespoke nurse call plate which can also fit the required number of data outlets. The design of the units ensures correct separation is maintained between electrical, data and medical gases.

The units can be fed either via a bespoke surface mounted vertical riser section, through the back of the unit through the main body or through the side of the unit through the end cap.

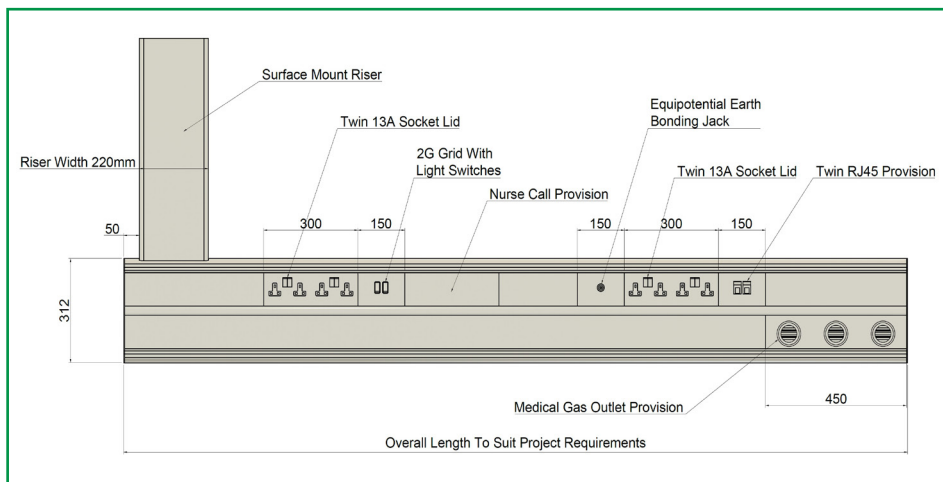

Specifications:

Mediline Horizontal Bedhead Unit

Electrical Socket Type	13A Switched/unswitched BS single or twin sockets - Quantity to client specification
Equipotential Sockets	POAG ID-6 as standard to client specification
Data Outlets	Single/Twin available to client specification
Gas Outlets	Quantities to suit project requirements - Pneumatech East/Zeus or Medaes gem 10 as standard
Nurse Call Provision	Compatible with all nurse call manufacturers
Construction	Powder coated aluminium extrusion with powder coated steel accessory plates to suit. Flexible to suit any number of beds with overall lengths designed to project requirements, supplied fully assembled with jointing kits where longer runs are required.
Service Entry	Back fed through carcass or side fed through end cap. A bespoke service riser is also available for surface service entry.
Colour	RAL9010 as standard. Optional colours available on request
Assembly	Supplied as fully assembled modules as standard with optional pre gassing and prewiring if required

T: 01244 288 299
enquiries@lipmedical.com

Office

Unit 60 AML Zone 2, Deeside Industrial Estate, Third Avenue, Deeside, Flintshire, CH5 2LA

