

# Mediwall Vertical Bedhead Unit


www.p-mgs.com



The future looks good.

# Vertical Bedhead Service Trunking


The Mediwall vertical bedhead unit offers a compact stylish solution which can be manufactured to suit each individual installation. Supplied as a fully assembled unit, electrical sockets, nurse call, data and medical gases are fitted to meet each individual client specification.

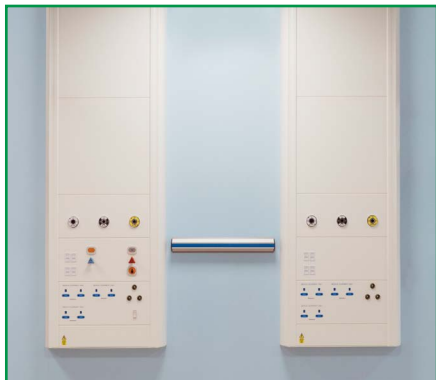

Electrical sockets are fitted and can be fully wired and tested with essential/non essential sockets available for standard ward areas and blue isolated earth sockets for A&E, critical care or operating departments. Medical gas terminal units from all manufacturers can either be provided as a fully pre-gassed and tested system or alternatively, mounting provision can be included for installation on site.

All types of nurse call can be easily integrated with provision for a bespoke nurse call plate which can also fit the required number of data outlets. The design of the units ensures correct separation is maintained between electrical, data and medical gases.

Accessories such as medical rail can be fitted with the option of high level sockets and data also available where patient monitoring is required.

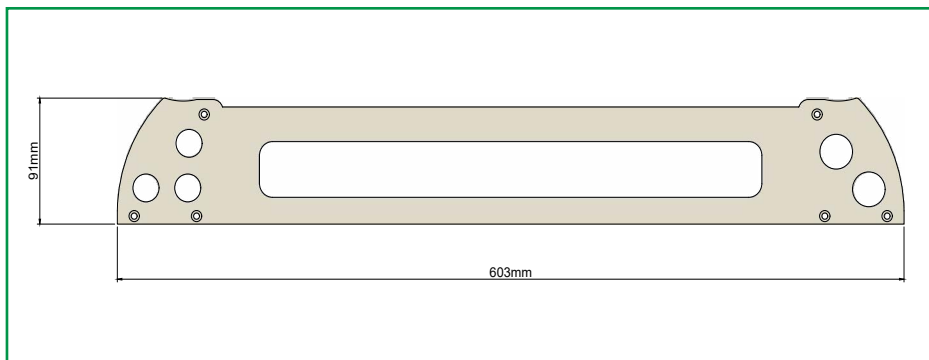

## Specifications:

# Vertical Bedhead Service Trunking

<b>Electrical Socket Type</b>	<b>13A Switched/unswitched BS single or twin sockets - Quantity to client specification</b>
<b>Equipotential Sockets</b>	<b>POAG ID-6 as standard to client specification</b>
<b>Data Outlets</b>	<b>Single/Twin available to client specification</b>
<b>Gas Outlets</b>	<b>Quantities to suit project requirements - Penlon East/Zeus or Medaes gem 10 as standard</b>
<b>Nurse Call Provision</b>	<b>Compatible with all nurse call manufacturers</b>
<b>Construction</b>	<b>Powder coated aluminium extrusion with powder coated steel accessory plates to suit. Length to suit project requirements. Services entry from the top of the unit with gland points for electrical/data services and separate medical gas entry provided</b>
<b>Colour</b>	<b>RAL9010 as standard. Optional colours available on request</b>
<b>Assembly</b>	<b>Supplied as fully assembled modules as standard with optional pre gassing and prewiring if required</b>

T: 01246 474 242  
martin.berry@p-mgs.com

### Office

Pneumatech Medical Gas Solutions, Telford Crescent, Staveley, Derbyshire, S43 3PF

