

'DIASYS' Minor Ops Light



www.lipmedical.com

The future looks good.



LIP MEDICAL
HEALTHCARE DESIGN




The DIASYS is a minor surgical light composed of a 4-LEDs module. There's no need for a switch, as soon as the DIASYS is moved, it will switch on automatically. At the end of the procedure, the light needs to be oriented towards the ceiling to be automatically switched off. Thanks to the patented Motion-Matic technology, the light will switch off automatically after 30 minutes without any movement.

The DIASYS becomes a resolutely modern light, oriented towards ecology and sustainable development and is available as either Wall, Ceiling or Mobile.


Ceiling Mounted

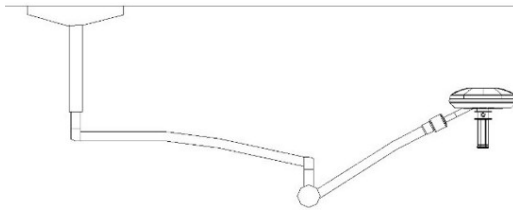

Wall Mounted

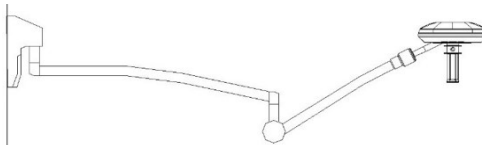

Mobile

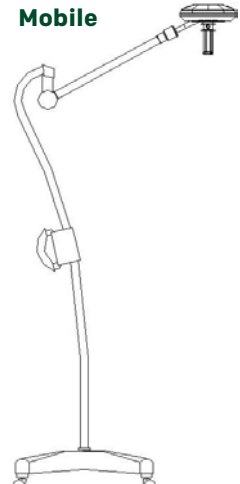

Specifications:

Minor Ops Lights

Light Source: LED

Maximum Central Lighting*: 65,000 LUX

Adjustable Luminous Intensity*: From 30,000 LUX to 65,000 LUX

Fixed Field Diameter*: 14cm

Depth of Field L1+L2 (60% Ec) 2008*: 128cm

Colour Rendering Index Ra*: 96

Red Colour Rendering R9*: 95

Colour Temperature*: 4200 K

LED Life Time: 60,000 h

Max. Consumption: 7w

Dimension of the Light: 30 x 24.5cm

T: 01244 288 299
enquiries@lipmedical.com

Office

Unit 60 AML Zone 2, Deeside Industrial Estate, Third Avenue, Deeside, Flintshire, CH5 2LA

